



# ADVANCE CARE *IP* Out-of-Bed Foldable Floor Mat

The Advance Care Out-of-Bed Foldable Floor Mat is a pressure sensitive device which sends off an alert the moment a resident steps out of their bed.

Composed of textured, non-slip material, the Out-of-Bed Foldable Floor Mat is designed to be both warm and comfortable underfoot. Aesthetically speaking, the mat is provided in a natural beige hue, ensuring that it is totally unobtrusive - simply becoming part of the room's decor.

With dimensions of 1200mm x 600mm, staff need not be concerned that a resident will exit their bed un-noticed; the length of the mat parallels with the length of a standard hospital bed.

The mat is water resistant, washable and durable, ensuring a lengthy product life.

A crucial addition in any facility's Dementia Care Management plan, when integrated with any Advance Care IP Nurse Call System, staff and management can specifically configure the Foldable Bed Exit Mat with a myriad of operation techniques. The Bed Exit Mat may also be configured to turn the lights on when the resident arises - a useful safety feature for those who are less mobile.


When the Foldable Floor Mat receives pressure, a high visibility reassurance light illuminates the 'CANCEL' button on the associated Call Button plate and the Call Button may also 'beep' periodically (if activated), providing reassurance to the user. The Out-of-Bed Foldable Floor Mat also has its own easily identifiable alert, such as 'Out of Bed Room 1'. This alert is also immediately displayed at the Nurse Station and sent to associated paging devices.

The addition of the Advance Care Out-of-Bed Foldable Floor Mat Interface allows for both the Floor Mat and a Call Pendant to be used simultaneously. The switch on the front of the Interface Device allows for the Foldable Floor Mat to be deactivated without needing to reposition the mat to avoid false alarms.

This product integrates with most Nurse Call Systems via a 6.35mm pendant plug and is a valued addition to any Advance Care IP Nurse Call installation.

**Features**

- Non slip
- Foldable
- Water resistant
- Pressure sensitive
- Simple plug-in use
- Durable construction
- Reliable out of bed detector
- Suitable with any level of care
  - Neutral in colour
  - Warm underfoot
  - Large size

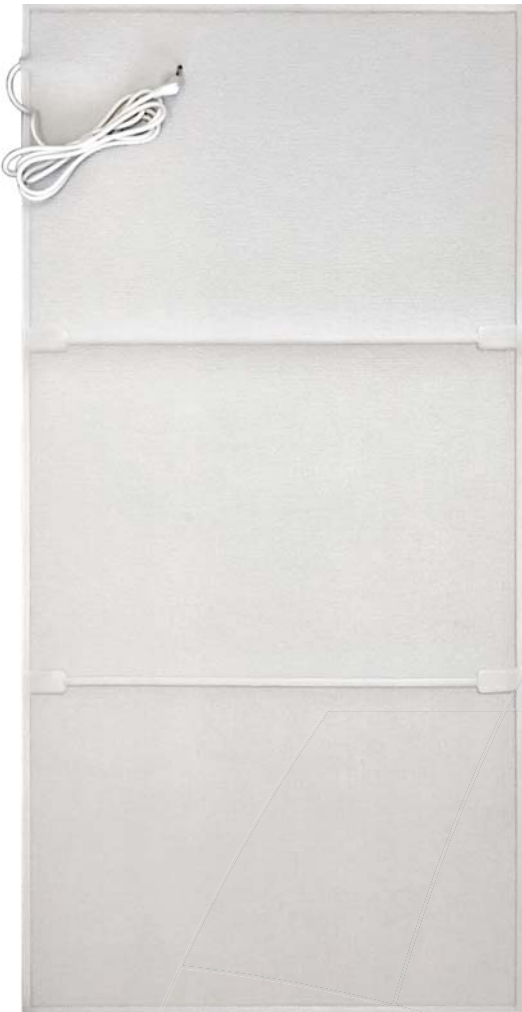

Code	Description
AC-FFM	Out-of-Bed Foldable Floor Mat
AC-FFM/I	Out-of-Bed Foldable Floor Mat Interface